

**BREAS**

UK Breathing Circuits and  
Consumables



## About Breas UK Ventilation Breathing Circuits

Breas UK assemble a range of breathing circuits for use with our ventilation device range. The circuits are assembled from components sourced by different suppliers, and are classified as “procedure packs” from a regulatory perspective.

In 2023, we made the decision to streamline our circuit portfolio, with the following aims in mind:

- To simplify things for our customers from both a procurement and usability perspective.
- To simplify things internally from a regulatory perspective by only using already CE marked components.
- To bring the circuit range up to date to be more suitable for use alongside our current ventilation device range.

Having taken 2 years to carefully phase many of our circuit variants, from 2026 onwards, Breas UK will only sell the configurations outlined in this brochure. For circuits intended for invasive use, or use with non-vented masks, we have made the decision to only use the Intersurgical purple leakage ports. These are widely used and well recognised around the UK, and hence we feel a safer option for patients.

These changes should make ordering the required circuits less confusing, and furthermore makes ongoing training around the use of Breas UK circuits simpler, both in the hospital and homecare settings. If needed, additional circuit components can be purchased separately when a circuit configuration needs adjusting by a clinician.

However, whenever a breathing circuit is changed, it is recommended that the patient’s ventilator parameters are reviewed and alarm settings are tested. This is to ensure alarm and therapy efficacy with the new circuit/ components (which may have different resistance, compliance and leakage properties to the previous ones used by the patient).

Please note, if a ventilator has a pre use test/learn circuit feature, it is always advisable to run this test whenever a circuit is changed, or as suggested to by the prescribing clinician.

As always, any questions, we are here to help.

BW  
The Breas UK Team


[salesUK@breas.com](mailto:salesUK@breas.com)



Page	Contents
1.	<u>Invasive 22mm Circuits - Adult (7 day use)</u>
2.	<u>Invasive 22mm Circuits - Adult (14 day use)</u>
3.	<u>Invasive 22mm Circuits - Paediatrics (7 day use)</u>
4.	<u>Invasive 15mm Circuits - Paediatrics (7 day use)</u>
5.	<u>Non-Invasive 22mm Circuits (7 day use)</u>
6.	<u>Non-Invasive 15mm Circuits (7 day use)</u>
7.	<u>22mm Non-Invasive Circuit (7 day use)</u>
8.	<u>Invasive 22mm Active Valve Circuits (7 day use)</u>
9.	<u>Miscellaneous Ventilation Circuits</u>
10.	<u>Ventilation Circuits - Additional Items</u>
11.	<u>Mouthpiece Ventilation Circuits/Kits</u>
12.	<u>15mm and 22mm CPAP tubing (6 month use)</u>
13.	<u>Nippy Clearway and Clearo circuits</u>
14.	<u>Clearway interfaces and consumables</u>


## Invasive 22mm Heated Leakage Circuit for use with MR850 Humidifier with Adult Catheter Mount (7 day use)

Product Code: 008056 | NHSSC Code: [FAG2545](#) | Box Quantity of 10

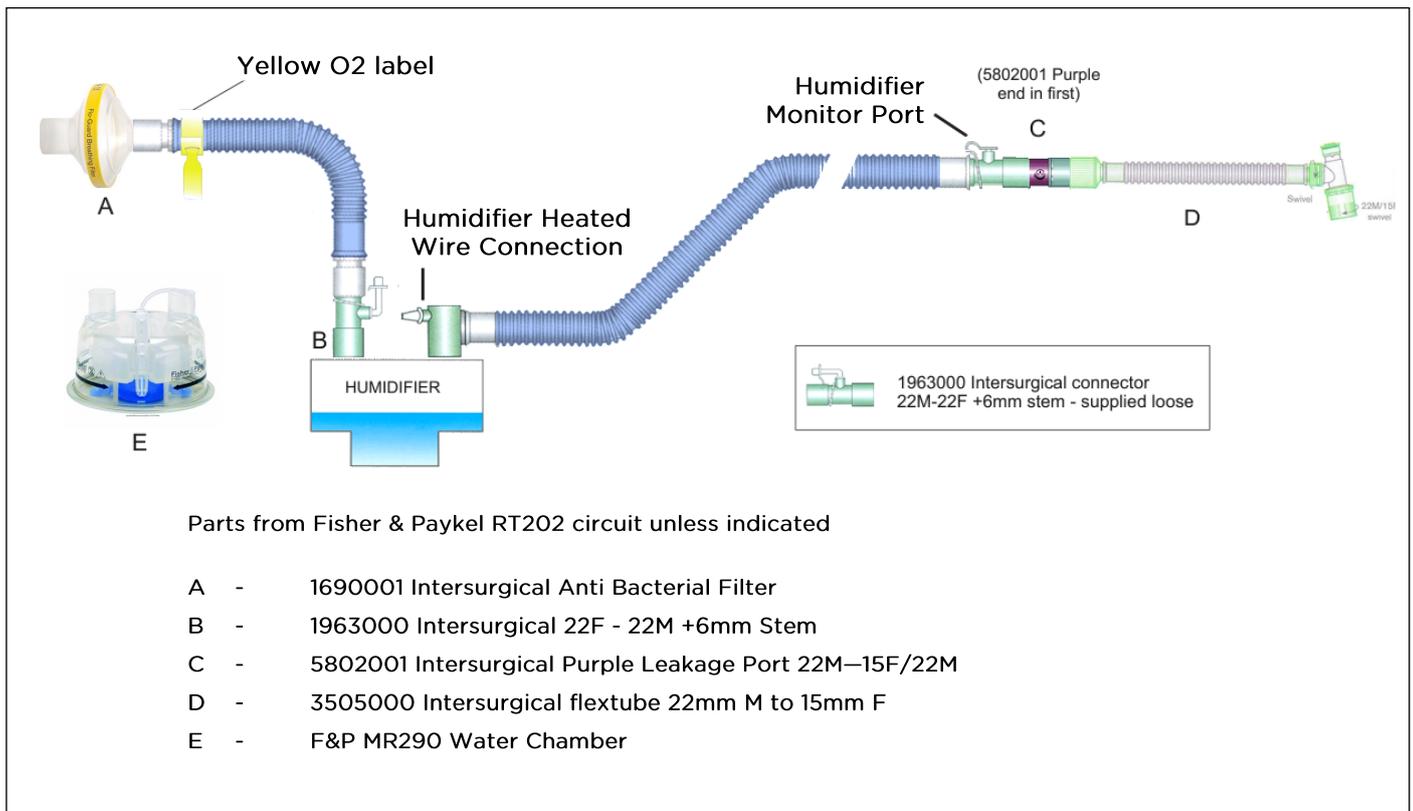

## Invasive 22mm Dry Leakage Circuit with Adult HME and Catheter Mount (7 day use)

Product Code: 008050 | NHSSC Code: [FAG2544](#) | Box Quantity of 10

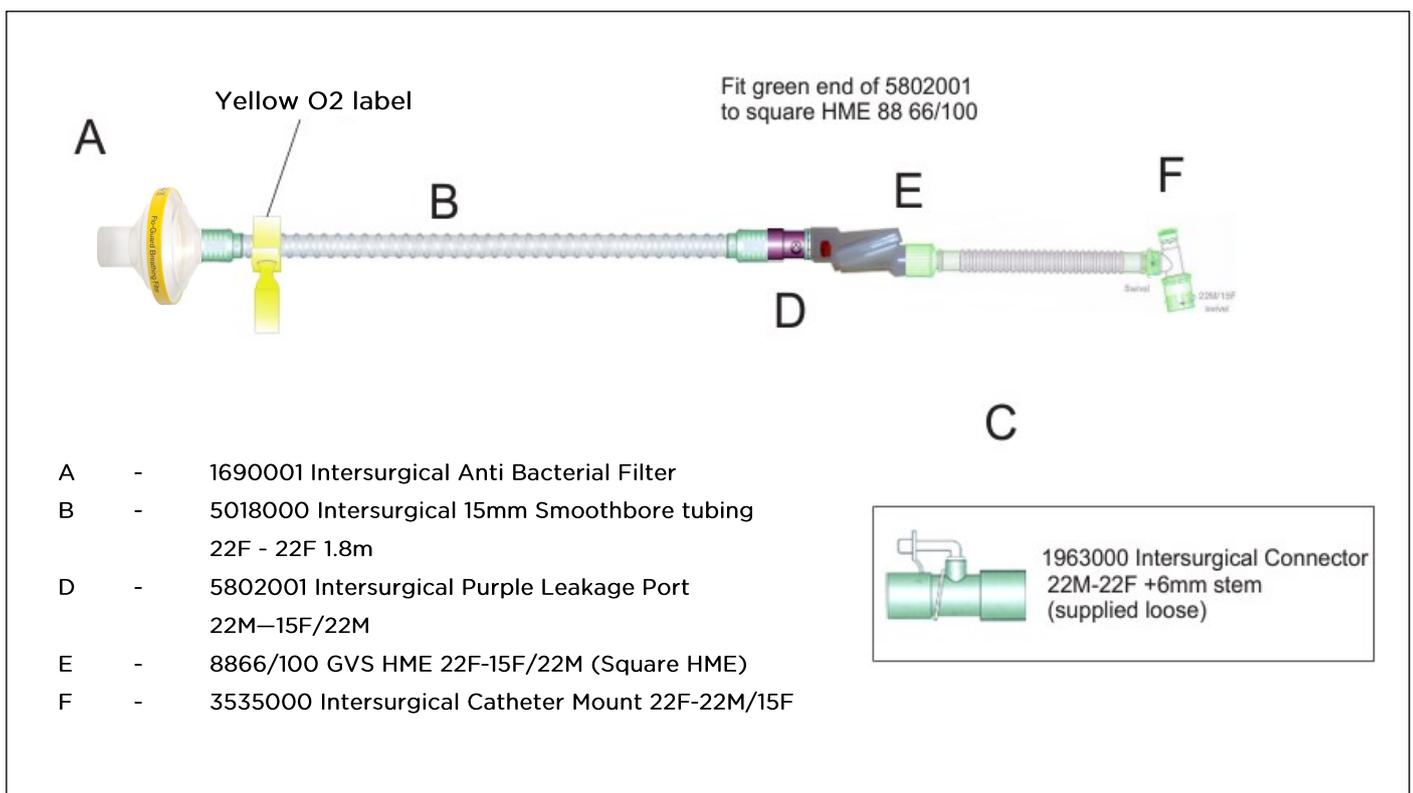

[BACK TO CONTENTS PAGE](#)


# Invasive 22mm Heated Circuit with Leakage Port and Adult Catheter Mount (14 day use) for use with External Humidifier

Product Code: 010471 | NHSSC Code: [FFH85100](#) | Box Quantity of 10

**Note:** This heated circuit is intended for 14 day usage before disposal, and is supplied with spare components which are to be changed after 7 days. This is one part of Breas UK's attempt to reduce overall carbon footprints for home ventilated patients in the wider drive towards net zero.


Spares supplied Loose (1 of each to be changed after 7 days)

<p>A - 1690001 Intersurgical Anti Bacterial Filter (7 day use)</p> <p>B - 51007091 Vincent Medical M 22mm Tubing, 0.6M length (14 day use)</p> <p>C - 51007091 - Vincent Medical Humidifier Temperature Monitor Point and Heated Wire Connection (14 day use)</p> <p>D - 51007091 - Vincent Medical 22mm Heated Tubing (14 day use)</p>	<p>E - 51007091 Vincent Medical Vent-X Membrane (14 day use)</p> <p>F - 51007091 Vincent Medical Humidity Monitoring Point (14 day use)</p> <p>G - 5802001 - Intersurgical Purple Leakage Port 22M—15F/22M (7 day use)</p> <p>H - 51002881 - Vincent Medical Auto Fill VHC20 Water Chamber (14 day use)</p>
---	---

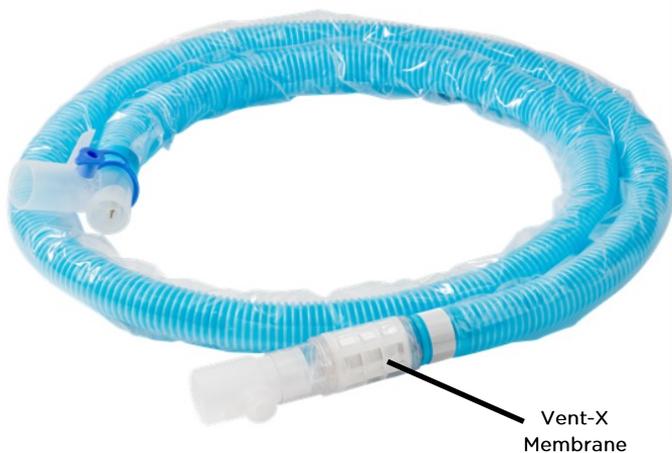

**This configuration from Breas uses a 22mm heated circuit which is intended for 14-day use.**


**The circuit has an additional polythene layer around the tubing, in order to regulate the ambient temperature around the heated part of the circuit. In turn, this controls unwanted condensation build up in the heated circuit.**

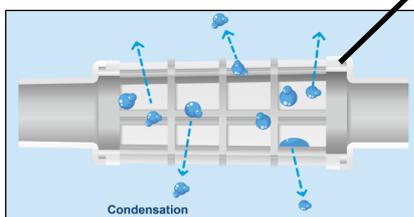

**It also possesses the “Vent-X” Membrane technology which allows excess moisture to diffuse out of the breathing circuit.**


## Invasive 22mm Humidified/Heated Leakage Circuit for use with MR850 Humidifier with Paediatric Catheter Mount (7 day use)

Product Code: 008059 | NHSSC Code: Not listed | Box Quantity of 10

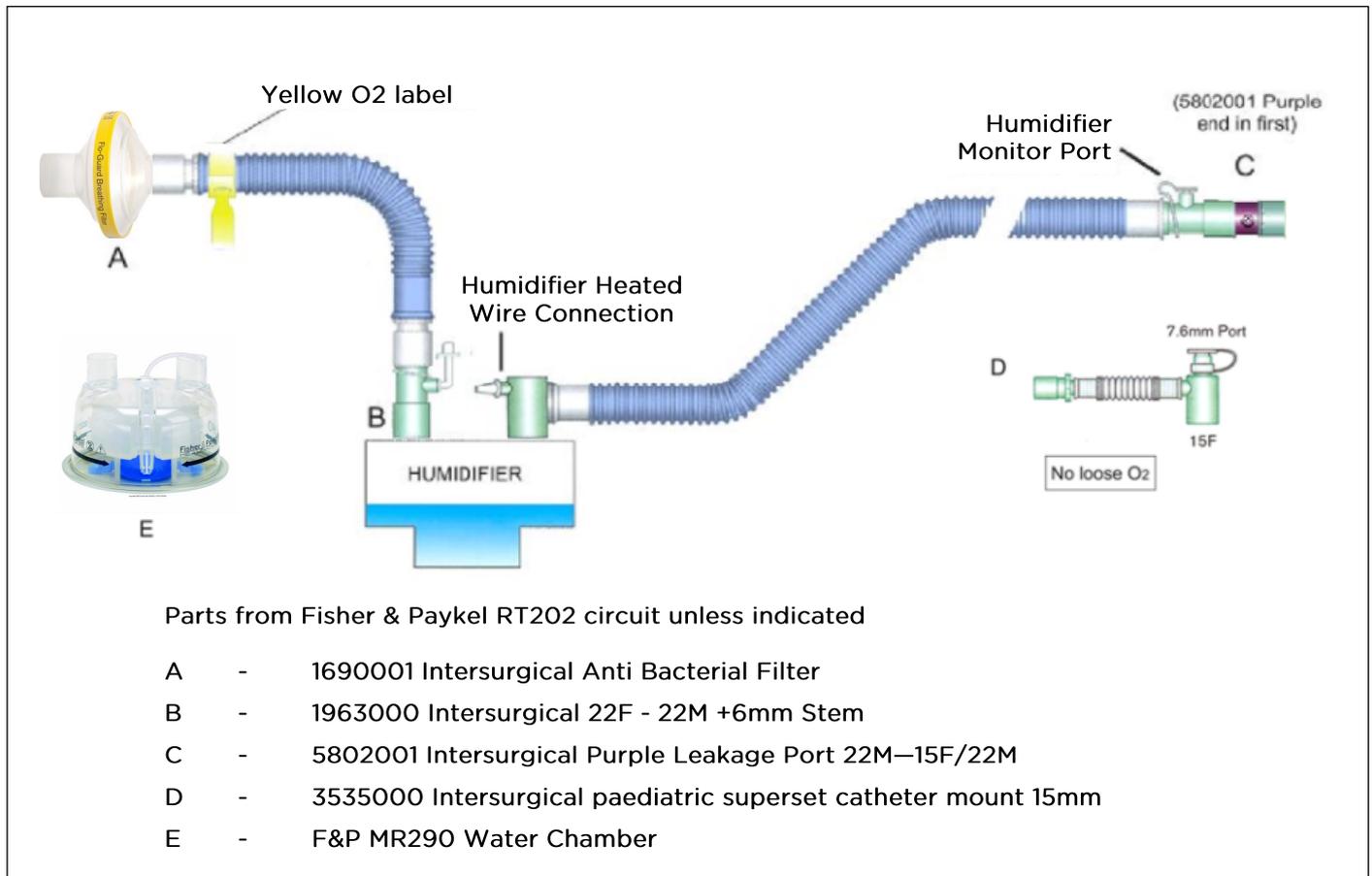

## Invasive 22mm Dry Leakage Circuit with Paediatric HME and Catheter Mount (7 day use)

Product Code: 008048 | NHSSC Code: Not listed | Box Quantity of 10

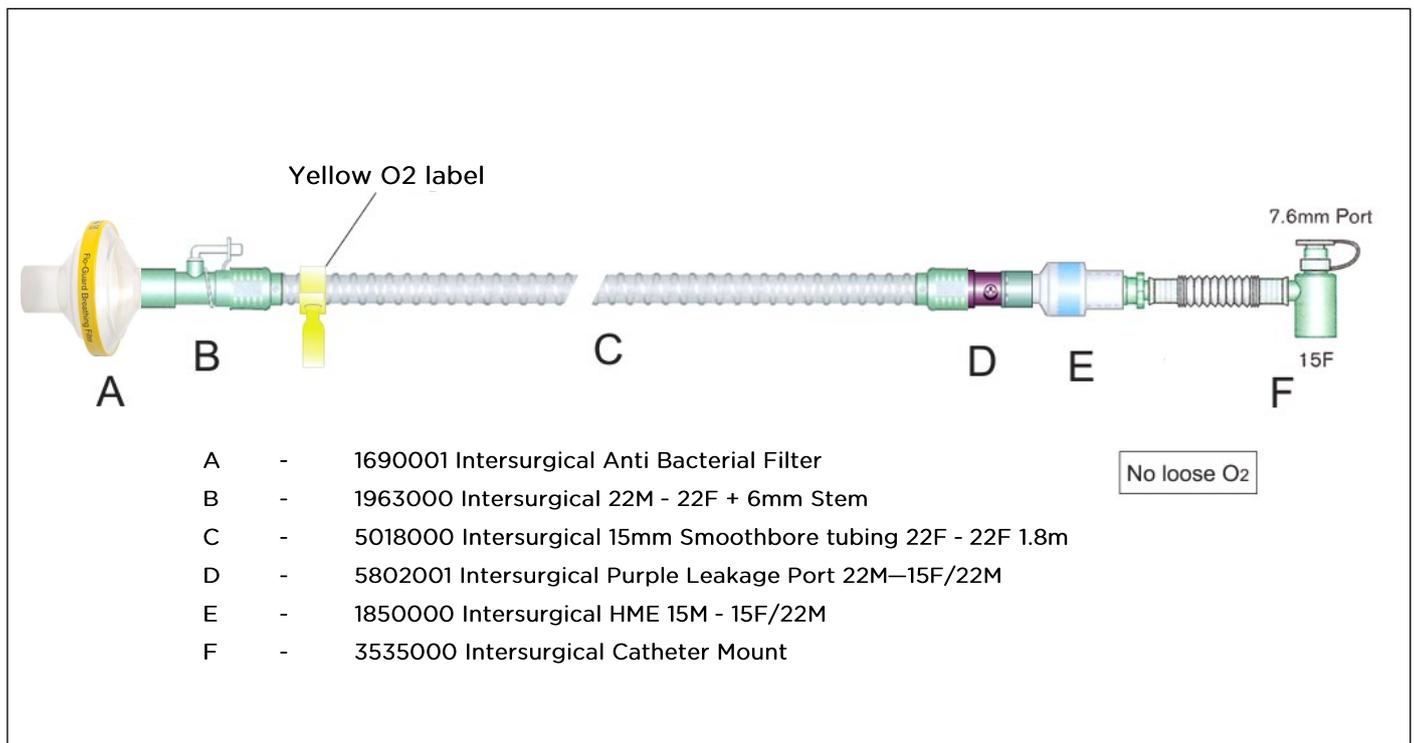

[BACK TO CONTENTS PAGE](#)

## Invasive 15mm Heated Leakage Circuit for use with MR850 Humidifier (7 day use)

Product Code: 009696 | NHSSC Code: [FFH85011](#) | Box Quantity of 10

Yellow O2 label

Humidifier Monitor Port

Humidifier Heated Wire Connection

INTERSURGICAL

A	-	1690001 Intersurgical Anti Bacterial Filter	E	-	Humidifier Temperature Monitor Point
B	-	5006000 Intersurgical 22mm Smoothbore Tubing (600mm length)	F	-	5802001 Intersurgical Purple Leakage Port 22M-15F/22M
C	-	Humidifier Temperature Monitor Point and Heated Wire Connection	G	-	Intersurgical Swivel Catheter Mount
D	-	5015008 Intersurgical 15mm Smoothbore heated Tubing (1500mm length)	H	-	2310000 Intersurgical Auto Fill Water Chamber

## Invasive 15mm Dry Leakage Circuit (7 day use)

Product Code: 008082 | NHSSC Code: [FFH85008](#) | Box Quantity of 10

Yellow O2 label

7.6mm Port

15F

A	-	1690001 Intersurgical Anti Bacterial Filter
B	-	5016000 Intersurgical 15mm Smoothbore tubing
C	-	1962000 Intersurgical 15F - 22F Connector
D	-	5802001 Intersurgical Purple Leakage Port 22M-15F/22M
E	-	1850000 Intersurgical HME 15M - 15F/22M
F	-	3535000 Intersurgical Catheter Mount

1963000 Intersurgical Connector 22M-22F +6mm stem (supplied loose)

## Non-Invasive 22mm Humidified Circuit for use with Vented Interface (7 day use)

Product Code: 008111 | NHSSC Code: [EFH85010](#) | Box Quantity of 10

Yellow O2 label

A B C D E F

HUMIDIFIER

No loose O2

F


- A - 1690001 Intersurgical Anti Bacterial Filter
- B - 1961000 Intersurgical Connector 22M-22F
- C - 1574000 Intersurgical 800mm (two sections)
- D - 1963000 Intersurgical 22F—22M + 6mm Stem
- E - 5018000 Intersurgical 22mm Smoothbore tubing
- F - F&P HC325 Manual Fill Water Chamber

## Non-Invasive 22mm Dry Circuit for use with Vented Interface (7 day use)

Product Code: 008024 | NHSSC Code: [FAG4538](#) | Box Quantity of 10

A B

Yellow O2 label

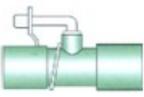

1963000 Intersurgical Connector  
22M-22F +6mm stem  
(supplied loose)

- A - 1690001 Intersurgical Anti Bacterial Filter
- B - 5018000 Intersurgical 22mm Smoothbore tubing

[BACK TO CONTENTS PAGE](#)

## Non-Invasive 15mm Humidified Circuit for use with Vented Interface (7 day use)

Product Code: 008113 | NHSSC Code: [FAG4321](#) | Box Quantity of 10

Yellow O2 label

A B C D E F

HUMIDIFIER

1962 is fitted in small end of 5016

No loose O2

F - F&P HC325 Manual Fill Water Chamber

A - 1690001 Intersurgical Anti Bacterial Filter

B - 5006000 Intersurgical 22mm Smoothbore tubing 0.6m

C - 1963000 Intersurgical 22F - 22M +6mm Stem

D - 5016000 Intersurgical 15mm Smoothbore tubing 1.6m

E - 1962000 22mmF to 15mmM connector fitted to small end tube

F - F&P HC325 Manual Fill Water Chamber

## Non-Invasive 15mm Dry Circuit for use with Vented Interface (7 day use)

Product Code: 008028 | NHSSC Code: [FAG4369](#) | Box Quantity of 10

1962 is fitted in small end of 5016

Yellow O2 label

A B C

1963000 Intersurgical Connector 22M-22F +6mm stem (supplied loose)

A - 1690001 Intersurgical Anti Bacterial Filter

B - 1962000 22mmF to 15mmM connector fitted to small end tube

C - 5016000 Intersurgical 15mm Smoothbore tubing 1.6m

## Non-Invasive 22mm Humidified/Heated Circuit for use with Vented Interface (7 day use)

Product Code: 008074 | NHSSC Code: [FFH85007](#) | Box Quantity of 10

Yellow O2 label

B

A

C

D

HUMIDIFIER

HUMIDIFIER Heated Wire Connection

Humidifier Monitor Port

E

1963000 Intersurgical Connector 22M-22F +6mm stem (supplied loose)

**WARNING: CIRCUIT DOES NOT HAVE EXHALATION PORT OR CATHETER MOUNT**

Parts from Fisher & Paykel RT202 circuit unless indicated

- A - 1690001 Intersurgical Anti Bacterial Filter
- B - 1961000 Intersurgical 22M-22F connector
- C - 1574000 Intersurgical 22mm tube 800mm = 2 SECTIONS
- D - 1963000 Intersurgical 22F - 22M +6mm Stem
- E - F&P MR290 Water Chamber

## Non-Invasive 22mm Dry Leakage Circuit for use with Non-Vented Interface (7 day use)

Product Code: 008026 | NHSSC Code: [FAG4442](#) | Box Quantity of 10

Yellow O2 label

B

A

C

Purple end in first

1963000 Intersurgical Connector 22M-22F +6mm stem (supplied loose)

- A - 1690001 Intersurgical Anti Bacterial Filter
- B - 5018000 Intersurgical Smoothbore Tubing 1.8M 22F-22F
- C - 5802001 Intersurgical Leakage Port 22M-22F

[BACK TO CONTENTS PAGE](#)

## Invasive 22mm Humidified/Heated IPPV Circuit for use with MR850 Humidifier (7 day use)

Product Code: 008101 | NHSSC Code: [FFH85009](#) | Box Quantity of 10

Parts from Fisher & Paykel RT202 circuit unless indicated

A	-	1690001 Intersurgical Anti Bacterial Filter
B	-	1963000 Intersurgical 22F - 22M +6mm Stem
C	-	Valve & Tubing (patient line) from Intersurgical 5194009
D	-	3505000 Intersurgical 22mm F to 15mm F
E	-	MR290 Water Chamber

## Invasive 22mm Dry IPPV Circuit (7 day use)

Product Code: 008034 | NHSSC Code: [FFH85004](#) | Box Quantity of 10

A	-	1690001 Intersurgical Anti Bacterial Filter
B	-	5194009 Intersurgical Valve, 22mm Tubing and pilot line
C	-	8866/100 GVS HME
D	-	3505000 Intersurgical Catheter Mount 22mm F to 15mm F



Pilot Line Connector for using Nippy 4+/Vivo 45LS with an active exhalation valve circuit

Product Code: 007654 | NHSSC Code: [FAG5679](#)

Supplied with every Nippy 4+, or sold individually

[BACK TO CONTENTS PAGE](#)

# Miscellaneous Ventilation Circuits



22mm SeptiFlow Patient Circuit with active exhalation valve

Product Code: 007616

NHSSC Code: [FAG4698](#)

Sold individually



Pilot Line Connector for using Nippy 4+ with an active exhalation valve circuit

Product Code: 007654

NHSSC Code: [FAG5679](#)

Sold individually



Heated Patient Circuit (Vivo 1-2-3)

Product Code: 006990

NHSSC Code: [FAG4730](#)

Sold individually



Heated Patient Circuit (Nippy 4 / 4+)

Product Code: 006193

NHSSC Code: [FAG4650](#)

Sold individually

[BACK TO CONTENTS PAGE](#)

# Ventilation Circuits - Additional Items



Anti Bacterial Filter

Code: 008018

NHSSC: [FFX85139](#)

Sold individually



Swivel Leakage Port

Code: 004426

NHSSC: [FAG4685](#)

Sold individually




HC325 Manual Fill Water Chamber

Code: 008351

NHSSC: [FDA85108](#)

Box Quantity of 10

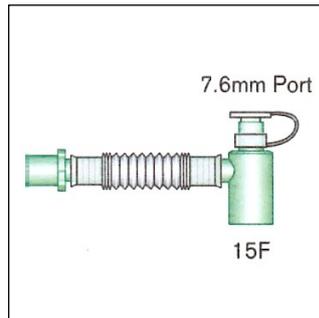

Extendable Paediatric Catheter Mount (with suction port)

Code: 008307

NHSSC: [FDI85044](#)

Sold individually




15mm Heated Moisture Exchange (HME)

Code: 008019

NHSSC: [FFX85124](#)

Sold individually

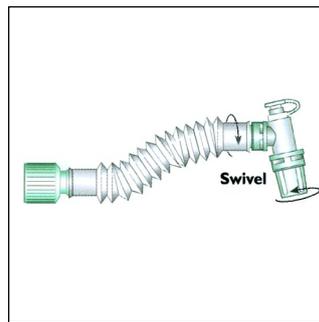

Extendable Swivel Catheter Mount (with suction port)

Code: 008306

NHSSC: [FDI85045](#)

Sold individually




22mm Heated Moisture Exchange (HME)

Code: 008021

NHSSC: [FFX85138](#)

Sold individually

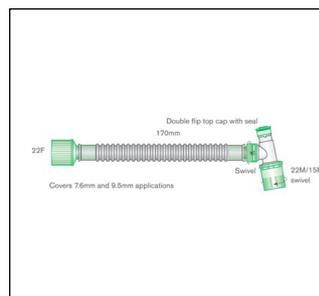

Extendable Swivel Catheter Mount (with suction port)

Code: 008302

NHSSC: [FDI85043](#)

Sold individually




Leakage Port

Code: 008130

NHSSC: [FFX85132](#)

Sold individually

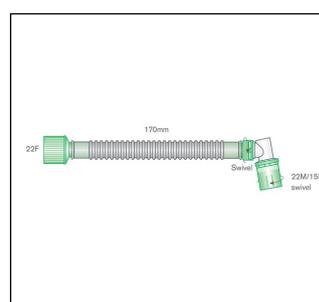

Extendable Swivel Catheter Mount (no suction port)

Code: 008303

NHSSC: [FDI85046](#)

Sold individually

# Mouthpiece Ventilation Kits - UK

The Nippy 4 can deliver an MPV-Pressure mode, and the Nippy 4+ and Vivo 65 ventilators are capable of delivering MPV-Pressure and Volume modes for your patients. You can find instructions in the User Manuals, and tutorial videos explaining the use of MPV can be found by clicking on or scanning the QR code. Breas have all of the required accessories available for patients using Mouthpiece Ventilation.



SCAN ME



# CPAP tubing



## 22mm Lightweight Smoothbore Breathing Circuit

Product Code: 008426

NHSSC Code: [FAG4908](#)

Not sold individually - only available in box quantity of 30

Discounts available for agreed annual volumes



## 15mm Lightweight Smoothbore Breathing Circuit

Product Code: 008427

NHSSC Code: [FAG4726](#)

Not sold individually - only available in box quantity of 30

Discounts available for agreed annual volumes

# Clearway Circuits




## MI-E circuit

Product Code: 008157 | NHSSC Code: [FAG2543](#) | Box Quantity of 10



## IPPB circuit

Product Code: 008185 | NHSSC Code: [FFX85123](#) | Box Quantity of 10

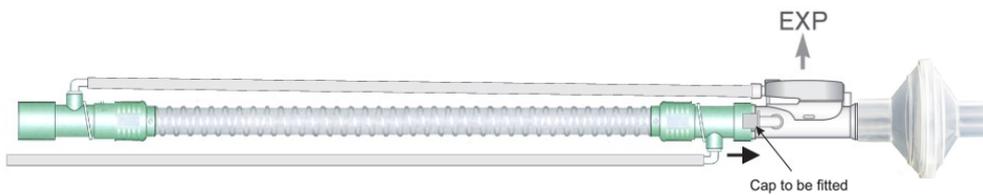

---

# Clearo Circuits



## MI-E circuit

Product Code: 008276 | NHSSC Code: [FAG6536](#) | Box Quantity of 10



## IPPB/NIV circuit

Product Code: 008277 | NHSSC Code: [FAG6530](#) | Box Quantity of 10



[BACK TO CONTENTS PAGE](#)

# Clearway and Clearo Interfaces and consumables



**Mask Size 6 (XL Adult)**

Product Code: 008342

NHSSC: [FAG5657](#)

Sold Individually



**Mask Size 5 (Large Adult)**

Product Code: 008341

NHSSC: [FAG5677](#)

Sold Individually



**Mask Size 4 (Adult)**

Product Code: 008340

NHSSC: [FAG5673](#)

Sold Individually



**Mask Size 3 (Small Adult)**

Product Code: 008339

NHSSC: [FAG5664](#)

Sold Individually



**Mask Size 2 (Paediatric)**

Product Code: 008338

NHSSC: [FAG5669](#)

Sold Individually



**Mask Size 1 (Infant)**

Product Code: 008337

NHSSC: [FAG5663](#)

Sold Individually



**Mask Size 0 (Neonate)**

Product Code: 008336

NHSSC: [FAG5668](#)

Sold Individually



**Mouthpiece**

Product Code: 008160

NHSSC: [FBU91349](#)

Box Quantity of 10



**Mouthpiece and Nasal Clip**

Product Code: 008161

NHSSC: [FBU91350](#)

Box Quantity of 10



**Catheter Mount**

Product Code: 008159

NHSSC: [FAG6528](#)

Box Quantity of 10



**Anti Bacterial Filter**

Product Code: 008169

NHSSC: [FAG6522](#)

Box Quantity of 10



**Oxygen Coupler**

Product Code: 008170

NHSSC: [FFX85133](#)

Box Quantity of 10



**CW2 In-let Filter, White**

Product Code: 009065

NHSSC: [FFX85137](#)

Pack of 5



**CW2 In-let Filter, Grey**

Product Code: 009066

NHSSC: [FFX85136](#)

Pack of 5



# Breas UK - Breathing Circuits Feedback Form



Medical suppliers need to capture market feedback on all of their products. If you use a Breas breathing circuit, please take a moment to complete this short questionnaire to provide us with your experiences and thoughts.

Just scan or click the above QR code to provide your feedback.

# Breas UK - Breathing Circuits

## Feedback Form



© 2026 Breas Medical Ltd - All rights reserved. Breas Medical Ltd reserves the right to make changes in specifications and features shown herein, or discontinue the product described at any time without notice or obligation. Contact your Breas representative for the most current information. Breas and the Breas logo are trademarks of Breas Medical AB.

Breas Medical Ltd. Unit A2 The Bridge Business Centre, Timothy's Bridge Road, Stratford upon Avon, Warwickshire. CV37 9HW. United Kingdom. +44(0)1789 293 460 Company registered in England and Wales No. 07190611.

Document Number STR-MAR-MAR-008019-Rev\_4